

50x100mm

SAFETY INSTRUCTIONS


Before using this massage device, please read the instruction carefully to ensure perfect function.

- Please avoid the equipment coming in contact with water, high temperatures as well as direct exposure to the sun.
- If you depend on electrical aids such as a pace maker, please consult your physician prior to usage in any case.
- By 2 x AA batteries (not included)
- Do not mix old and new batteries
- Do not dispose of batteries in fire, batteries may explode or leak
- Do no charge non-chargeable batteries. Remove batteries from product if not in use for an extended period of time

Operating instructions:

- Carefully remove the item from the packaging
- Open the battery cover & insert 2 x AA batteries. Being careful to observe proper polarity. Fully tighten.
- This unit has two speeds:
Press the button once to "Low" - lower the speed
Press the button again to "High" - maximum speed
Press the button again to "Off"

IMPORTED BY:
S & S ENTERPRISES LTD
26 PERIVALE PARK,
PERIVALE, UB6 7RL

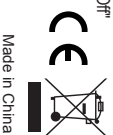

Made in China