

Rechargeable Neck Therapy Massager

H603 v01: 03/03/21

Please read and retain these instructions for future reference

vivadia



Dimensions

(H)5.8cm x(W) 17.5 x(D)14cm

Composition

ABS

1400mah battery (included)


Product Feature:

- Ergonomic design for a comfortable fit
- Wireless controller with large LCD display.
- Heating with low and high temperature setting.
- Low frequency acupuncture, 16 levels intensity and 6 types impulses.
- Vibration massage, 6 levels intensity, 6 kind of massage ways.
- Magnetic therapy, to improve micro circulation.
- Time setting from 5 min to 30 min adjustable.
- USB charged (cable included)
- One pair of pads connected to massager to send impulse to any part of the body.

ASSEMBLY INSTRUCTIONS / INSTRUCTIONS FOR USE:

Unpack the components carefully. Please check all parts have been removed from the packaging.


Parts

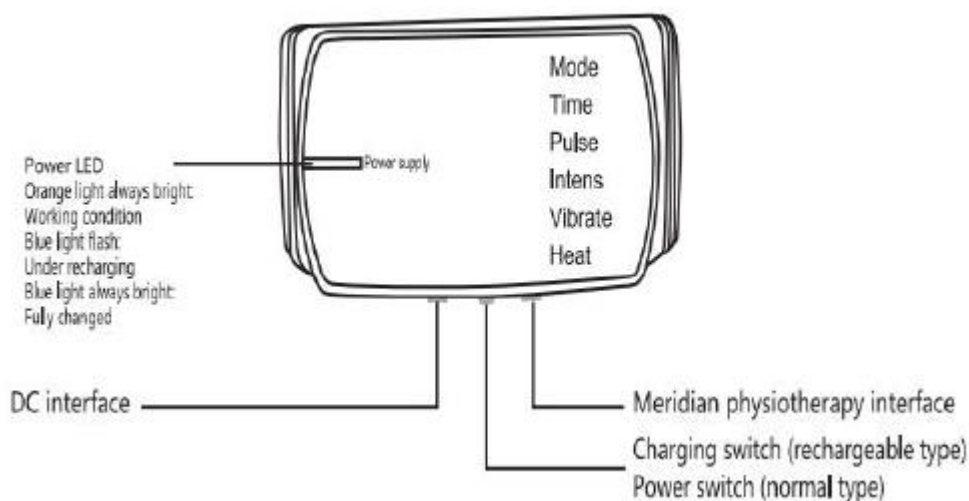


Rechargeable Neck Therapy Massager H603 Instructions for Use

Please read and retain these instructions for future reference

v001: 10/02/21

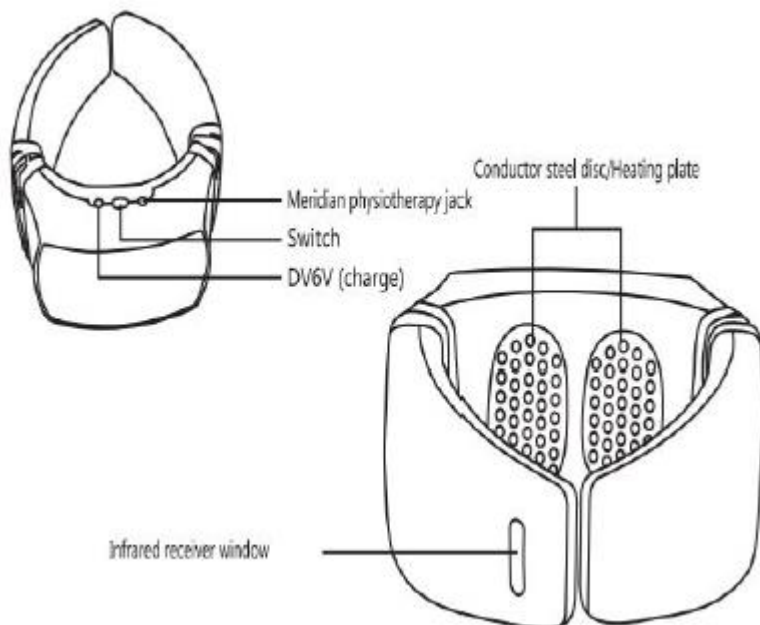

Accessories

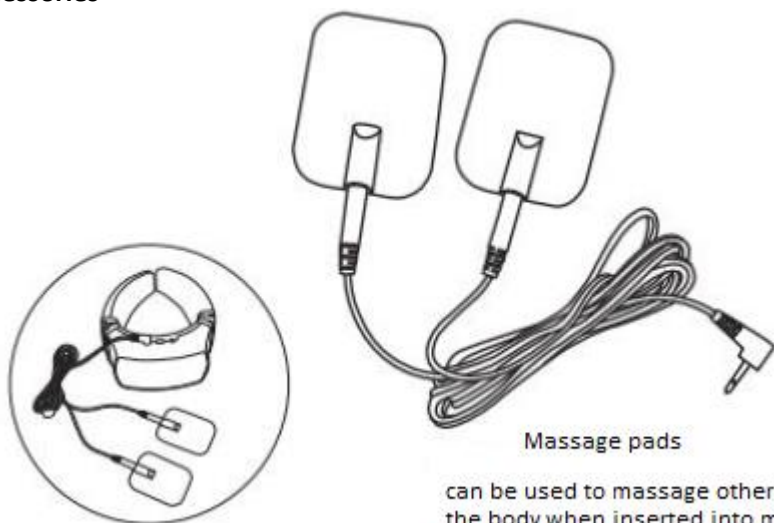


*Calls cost 7 pence per minute plus your telephone company's network access charge

The Enterprise Department
Head Office: 11 Bridge Street
Bishop's Stortford, CM23 2JU
Tel: 0844 482 4400*
Fax: 01279 756 595
www.coopersofstortford.co.uk

Rechargeable Neck Therapy Massager H603 Instructions for Use


Please read and retain these instructions for future reference

v001: 10/02/21

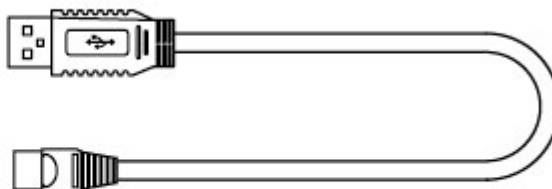

Power adapter
Insert into DC connector for charging

Remote control

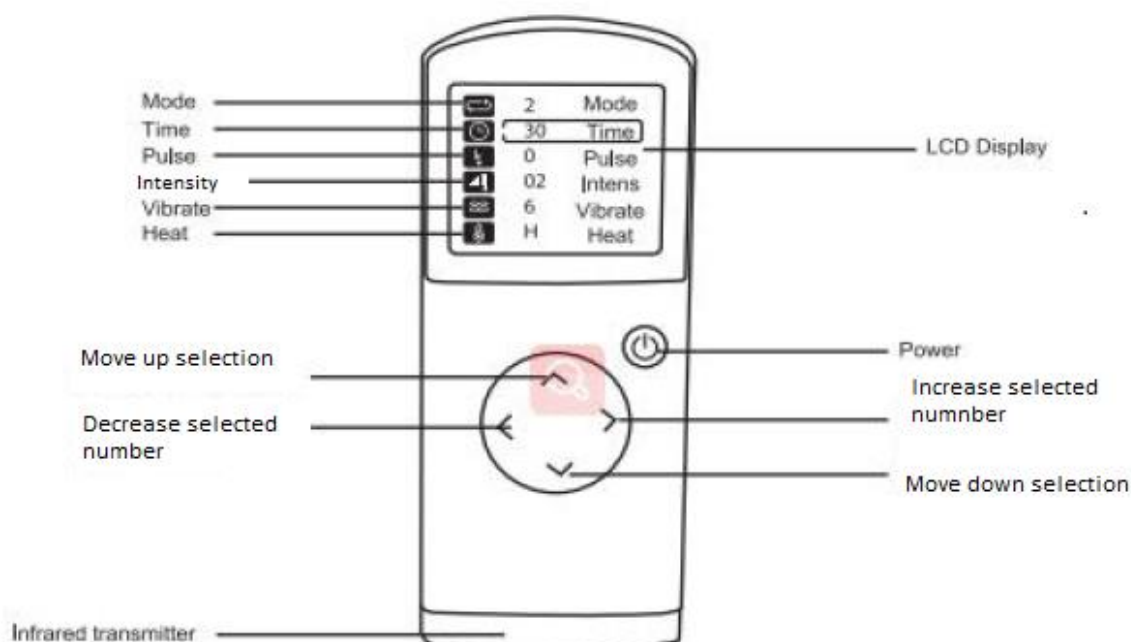

Method of Operation

Power Supply

Please plug the USB cable into a USB port to charge.

2. Built-in Rechargeable Battery Power Supply

The rechargeable battery supplies power without the massager being plugged into a USB port.

*Calls cost 7 pence per minute plus your telephone company's network access charge

The Enterprise Department
Head Office: 11 Bridge Street
Bishop's Stortford, CM23 2JU
Tel: 0844 482 4400*
Fax: 01279 756 595
www.coopersofstortford.co.uk

Rechargeable Neck Therapy Massager H603 Instructions for Use

Please read and retain these instructions for future reference

v001: 10/02/21

Function Operation

1. Firstly put the power switch to "ON" position. Secondly wet the conductor steel disc and your neck with a damp towel.
2. Thirdly put the massager close to the neck. Please use on clean skin.
3. Push the power button "begin" and the instrument buzzer will make a "click" sound. The power indicator lights up and the massager starts to work. The remote control LCD is turned on to enter auto working mode.

Set to


A:Heating Therapy

B:Full-featured Vibration Massage

C:First Level Vibration Intensity

D:Auto Timer 15 Minutes

Suggested use

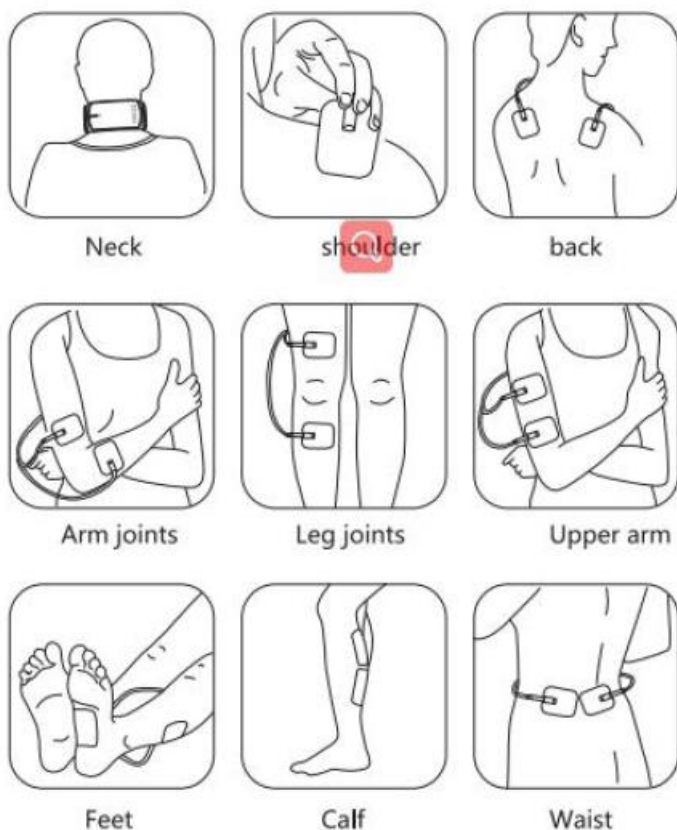

*Calls cost 7 pence per minute plus your telephone company's network access charge

The Enterprise Department
Head Office: 11 Bridge Street
Bishop's Stortford, CM23 2JU
Tel: 0844 482 4400*
Fax: 01279 756 595
www.coopersofstortford.co.uk

Rechargeable Neck Therapy Massager H603 Instructions for Use

Please read and retain these instructions for future reference

v001: 10/02/21


CARE AND CLEANING INSTRUCTIONS:

1. If the instrument stained with dirt and stains, please wipe with a soft cloth or dry towel.
2. After each use, please clean the electrode slice with a soft cloth or wrung-out damp towel.
3. If the electrode slice is very dirty, please use a soft cloth dampened with 75% medicinal alcohol.
4. The product should be stored in a cool dry place avoiding direct sun light.
5. No user serviceable parts.
6. Keep out of reach of children.

DISPOSAL

- Coopers of Stortford use recyclable or recycled packaging where possible.
- Please dispose of all packaging, paper, cartons, packaging in accordance with your local recycling regulations.
- Plastics, polybags – this is made from the following recyclable plastic.



Code & Symbol	 OTHER
Type of Plastic	ABS
Commonly used for	CD's, crisp packets, various flexible packaging, baby bottles, sunglass lenses
Notes	Not often or widely recycled

- At the end of the product's lifespan please check with your local council authorised household waste recycling centre for disposal.



vivadia

*Calls cost 7 pence per minute plus your telephone company's network access charge

The Enterprise Department
Head Office: 11 Bridge Street
Bishop's Stortford, CM23 2JU
Tel: 0844 482 4400*
Fax: 01279 756 595
www.coopersofstortford.co.uk