



Eye Massage Wand




Instructions for use:

1. Before use, charge the product using the USB cable included via the USB port as indicated in the diagram below.
2. Once charged, apply your usual face or eye cream to the upper and lower eyelid.
3. Using the on / off / mode button as indicated in the diagram below, press and hold to turn the product on. The product will light up red to indicate the heated setting, to change to the non-heated setting simply press the button again.
4. Select the desired mode and massage around the upper and lower lid ensuring to continuously move the product around.
5. When finished, hold down the on / off / mode button to turn the product off.
6. Using a damp towel clean the massage head and dry.

Specifications:

Power	$\leq 1W$
Vibration Frequency	12000rpm
Operating Current	20-500mA
Operating Voltage	3.7V
Charging Voltage	5V
Charging Time	Approx. 2 hours
Usage Time	Approx. 2 hours

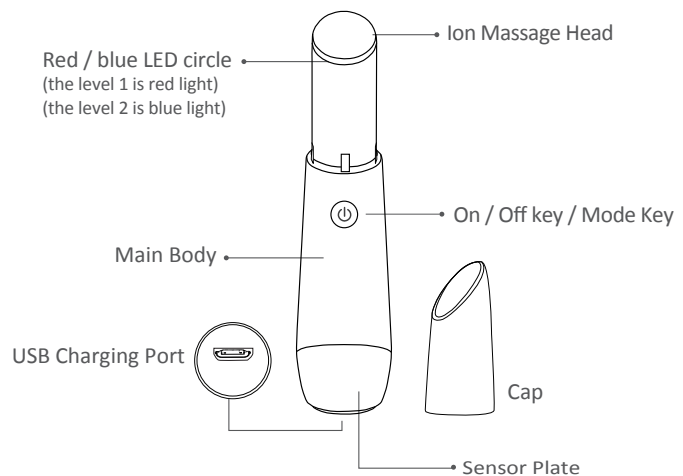

Caution:

- Test this product on a different area of the body such as the arm before use to check for any sensitivity or reaction.
- Do not wash this product in water, allow to cool after use and wipe clean with a damp tissue before drying.
- Always replace the cover after use.
- Store this product in a dry, ventilated area out of direct sunlight.
- When using this product, it is important to keep it moving and not hold the product in one place.
- If there is a suspected fault with the product, we do not recommend disassembling the product yourself, please consult a qualified technician or return the product to where you made your purchase from.

Warning:

- Always consult a doctor before use if you have pre-existing illnesses, have a pacemaker fitted or are pregnant.
- This product is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities or lack of experience or knowledge unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Never immerse the product in water.
- Never use the product outside.
- The product should be unplugged when not charging.
- Never try to turn on or use the product whilst it is charging.
- Do not plug in or use the product if the charging cable is damaged.

Please retain this instruction manual for future reference.



Product Code: N222342
Made in China for S&S Enterprises Ltd. 26 Perivale Park Middlesex UB7RL

