



## 3 in 1 Handheld Massager



Thank you for purchasing our 3 in 1 Handheld Massager.

We hope you are happy with your purchase.

**Please read the instructions carefully before use.**

### **Each set includes:**

- 1 x Massager Device
- 3 x Interchangeable Attachments (Copper Head, Cone Head, Fork Head)
- 1 x USB cable

### **Instructions for use:**

1. Select the appropriate attachment you wish to use and gently screw the attachment on to the massage gun head.
2. Press and hold the power button located on the back of the device to turn on.
3. Massage the required body parts at the desired speed and apply pressure as required.
4. After use, allow the device to cool down before gently unscrewing the attachment.


Whilst using ensure you do not over flex the unit as excessive bending may damage the device.



**Vertical Use**



**No Forward Bending**

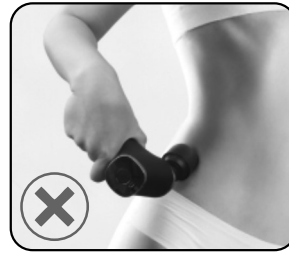

**No Upward Bending**



**No Backward Bending**

### **Charging:**

- This is a rechargeable unit for wireless use.
- Approximate charging time: 90 minutes.
- To charge simply plug in the USB charging cable to the charging outlet located next to the power button.
- **Do not** use the device whilst it is charging.

### **Maintenance:**

1. To clean and maintain the product, use a damp cloth to clean the device and dry completely before use.
2. Do not use the product for more than 30 minutes continuously, this may cause the device to overheat. If this happens, switch the device off and do not use or charge for 30 minutes to allow it to cool down.

### **Technical Specifications:**

<b>Dimensions</b>	16.4x15x4.8cm
<b>Material</b>	ABS+PVC
<b>Power supply</b>	Lithium Battery
<b>Output voltage/current</b>	Max 5.5W
<b>Battery capacity</b>	18650 1200 mAH
<b>Working time</b>	30 mins
<b>Operation</b>	one-button operation
<b>Speed</b>	2 speed user-friendly
<b>Db</b>	60dB
<b>Rotate speed</b>	1st 5300rpm 2nd 5700rpm

## **Warnings**

- Always consult a doctor before use if you have pre-existing illnesses, have a pacemaker fitted or are pregnant.
- This product is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities or lack of experience or knowledge unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children should be supervised at all times whilst using the product.
- Never immerse the product in water.
- Never use the product outside.
- The product should be unplugged when not charging.
- Never try to turn on or use the product whilst it is charging.
- Always give the product time to cool before changing or removing attachments.
- Never carry the product by the charging cable.
- Do not plug in or use the product if the charging cable is damaged.



Product Code: H22334

Made in China for S&S Enterprises Ltd, 26 Perivale Park Middlesex UB7RL

