




Rechargeable Handheld Massager with 6 Vibration Modes



Thank you for purchasing our Rechargeable Handheld Massager with 6 Vibration Modes.

We hope you are happy with your purchase.

Please read the instructions carefully before use.

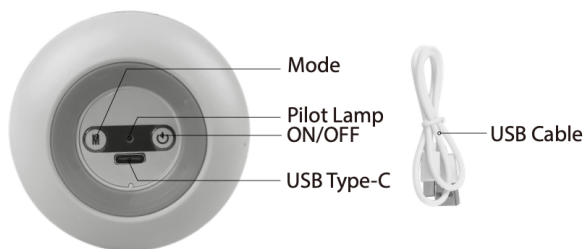

Each set includes:

- 1 x Massager Device
- 1 x USB cable

Instructions for use:

1. Press and hold the power button located on the back of the device to turn on (please see image above)
2. Massage the required body parts at the desired speed and apply pressure as required. The massager has 6 speeds, to switch between them simply press the power button once.
3. To turn the device off, hold the power button down.

Charging:

- This is a rechargeable unit for wireless use.
- Approximate charging time: 80 minutes.
- To charge simply plug in the USB charging cable to the charging outlet located next to the power button.
- When charging process begins, the light flashes blue, when fully charged, the light goes off
- **Do not** use the device whilst it is charging.

Maintenance:

1. To clean and maintain the product, use a damp cloth to clean the device and dry completely before use.
2. When it is not in use, you should store the product in a dry, clean location which is out of reach of children and animals
3. During a lengthy period of storage, check the charge status of the battery regularly and re-charge the battery if necessary in order to prevent a deep discharge.

Technical Specifications:

L7.2xcm x W7.2cm x D11.2cm

Code: H232522

DC IN: 5V-1A, 350mAh

Warnings

- Always consult a doctor before use if you have pre-existing illnesses, have a pacemaker fitted or are pregnant.
- This product is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities or lack of experience or knowledge unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Never immerse the product in water.
- Never use the product outside.
- The product should be unplugged when not charging.
- Never try to turn on or use the product whilst it is charging.
- Always give the product time to cool before changing or removing attachments.
- Never carry the product by the charging cable.
- Do not plug in or use the product if the charging cable is damaged.

Please retain this instruction manual for future reference.

Product Code: H232522

Made in China for S&S Enterprises Ltd. 26 Perivale Park Middlesex UB7RL

