3 IN 1 ACUPRESSURE PEN

Item No: M730 / 626670

Instructions for Use

v001: 13/05/2024

Java
wellness

ENGLISH INSTRUCTIONS

Power Details: Battery Type & details: 1xAA
Voltage details: DC 1.5V
Included / Not Included: Not Included


**ALWAYS USE ALKALINE BATTERIES.
DO NOT USE LITHIUM BATTERIES.**

Accessories: N/A

Please read this instruction manual thoroughly before starting and ensure you are familiar with the operation of your new 3 IN 1 ACUPRESSURE PEN

The meridian pen adopts human bioelectric pulse therapy and automatic acupoint seeking massage. The electric pulse stimulates the acupoints to stimulate cells and dredge the meridians. The current is conducted by hand and forms a backflow with the front section of the pen head of the meridians. Therefore, it is necessary to contact the backflow sheet at the handle of the machine. There is a slight numbness in the position of the hand. This is the normal phenomenon of bioelectric backflow of human body. Increasing the area of holding the reflux sheet can appropriately reduce the numbness.

[correctly use the acupuncture pen for others] while the helper holds the acupuncture pen, the other hand should contact the body skin of the user, otherwise the bioelectric circuit cannot be formed: The current intensity increases step by step from low to high (digital display 0 - no intensity, standby for 3-5 seconds and automatic shutdown; display 1 – lowest intensity, display 9 - maximum intensity) until it reaches the intensity acceptable to the user.

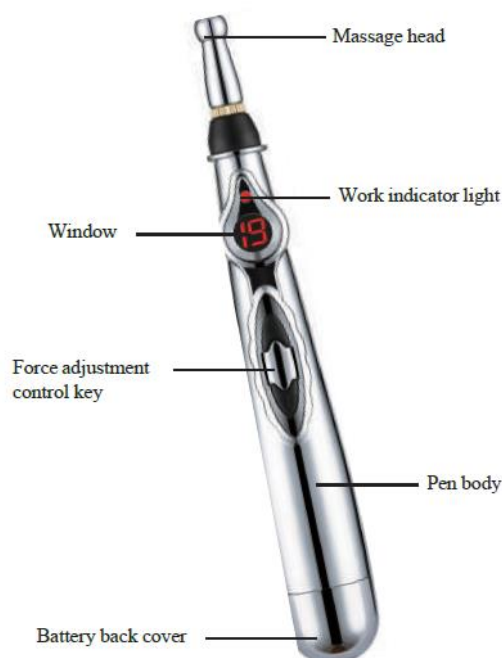1. Usage.

- 1) First, open the battery cover (as shown in the figure) and install the 1.5v 1xAA battery into the battery slot, and then tighten the back cover.
- 2) Use on skin:

- a) Before use, apply the gel to the body area, directly reaching the deep layer of the skin, locking in moisture, lubricating the skin, and allowing the three metal heads of the round cap, ball, and point to slide freely on the skin, increasing comfort.
- b) When it is turned on, press the "+" key, and the display will show "0". When it is pressed to 1, the indicator is flashing, indicating that it is working, then hold the machine reflux sheet, gently touch the probe of the pen on the human skin, and slowly slide and massage.
- c) The number has 0 - 9 ten numbers to show the strength of electric pulse scraping ("1" is the weakest and "9" is the strongest ") to increase the size, just press " + " " - "
- d) After use, please press the "0" key back and it will go out in 5 seconds.

3) Laser usage

When on, press and hold the "+" key for 3 seconds. At this time, the laser flashes, indicating that the laser is working Press and hold the "" - "key for 2 seconds. At this time, the laser is off

2. Precautions

The safety category of this product is BF type of internal power supply

- 1) Environmental protection: for relevant wastes (such as failed batteries and scrapped hosts), please don't throw them everywhere.
 - 2) When not in use for a long time, take out the battery to prevent liquid leakage from damaging the host
 - 3) Please do not open the host to repair it by yourself, otherwise the company will not be responsible for any consequences
 - 4) Never use the following electronic medical devices
 - (1) Pulse regulator and other in vivo transplanted medical electronic instruments
 - (2) Artificial heart, lung and other life sustaining medical electronic devices
 - (3) ECG recorder and other wearable medical electronic instruments.
 - 5) Under the following circumstances, please consult your doctor before use
 - (1) Patients with acute diseases
 - (2) Patients with malignant tumors
 - (3) Patients with infectious diseases
 - (4) Pregnant women
 - (5) Patients with severe heart disease
 - (6) Patients with abnormal blood pressure
 - (7) Patients with skin perception disorder or abnormal skin.
 - (8) Other patients who are receiving treatment and have caused accidents, especially when they feel physically abnormal.
- Be careful
- (1) Do not use with electronic devices.
 - (2) Do not use it for children with unconscious or unconscious expression. There is a risk of accidents or physical discomfort
 - (3) Do not use it in the bathroom or other places with high humidity, or use it while bathing, it will be strongly stimulated
 - (4) Do not use while driving. There is a risk of accidents and trouble caused by strong stimulation.
 - (5) Do not use if it irritates the skin and targets other than meridians

3. Main technical parameters

- 1) Power supply: DCL 5V; One No. 5 battery
- 2) Output: 3.7v; 300mA $\pm$ 50mA
- 3) Pulse frequency: 0.01 -300hz
- 4) Pulse width: 100ps-320 PS
- 5) Open circuit output electronic pulse voltage peak: 0-360V
- 6) Maximum output amplitude: 500 Ω load voltage is 3.3V

4. Equipment composition

Main engine, metal head and fuselage.

5. Instrument cleaning, disinfection and sterilization

1) Instrument cleaning

Remove the battery.

Use a soft brush to remove dust, hair and debris.

Wipe the surface with a clean, damp cloth that does not drip. It is recommended that the instrument be thoroughly cleaned after multiple use.

2) Disinfection and sterilization

Wipe the surface of the instrument repeatedly with cotton wool soaked in medical alcohol. After wiping, let it dry until the alcohol smell dissipates.

WARNINGS AND SAFETY INSTRUCTIONS:

WHEN USING ANY APPLIANCE, BASIC SAFETY PRECAUTIONS SHOULD ALWAYS BE FOLLOWED:

- READ ALL INSTRUCTIONS BEFORE USE.
- FOR DOMESTIC INDOOR.
- NOT FOR COMMERCIAL USE.
- BEFORE USE, CHECK THOROUGHLY FOR ANY DEFECTS AND DO NOT USE IF DEFECTS ARE FOUND. TAKE CARE NOT TO DROP THE APPLIANCE AS HEAVY IMPACTS MAY CAUSE INTERNAL DAMAGE.
- THE USE OF ANY ACCESSORY OR ATTACHMENTS NOT RECOMMENDED BY COOPERS OF STORTFORD MAY CAUSE INJURIES AND INVALIDATE ANY WARRANTY YOU MAY HAVE.
- DO NOT PLACE ON OR NEAR HEAT SOURCES.
- DO NOT COVER THE ITEM IN ANY WAY WHEN IN USE AS THIS MAY CAUSE OVERHEATING.
- DO NOT USE THIS APPLIANCE FOR ANY OTHER USE THAN THE INTENDED USE IT IS DESIGNED FOR.
- DO NOT LEAVE SWITCHED ON OR UNATTENDED.

BATTERIES

- **ALWAYS USE ALKALINE BATTERIES. DO NOT USE LITHIUM BATTERIES.**

- **BATTERY WARNING: KEEP OUT OF REACH OF CHILDREN. IF SWALLOWED OR PLACED INSIDE ANY PART OF THE BODY SEEK IMMEDIATE MEDICAL ATTENTION.**
- INSTALL BATTERIES IN ACCORDANCE WITH ITS POLARITY SYMBOLS. DO NOT FIT BATTERIES IN REVERSE.
- DO NOT SHORT CIRCUIT THE BATTERY TERMINALS.
- REMOVE THE BATTERIES IF LEFT FOR LONG PERIODS.
- **KEEP BATTERIES AWAY FROM CHILDREN AND PETS.**
- DO NOT MIX USED AND NEW BATTERIES.
- **NEVER DISPOSE OF BATTERIES IN A FIRE.**
- BATTERY DISPOSAL, SPENT OR EXPIRED BATTERIES MUST BE DISPOSED OF AND RECYCLED IN COMPLIANCE WITH YOUR LOCAL AUTHORITY REGULATIONS.
- FOLLOW THE BATTERY MANUFACTURER'S SAFETY, USAGE, AND DISPOSAL INSTRUCTIONS.

CARE AND CLEANING INSTRUCTIONS

- ALWAYS SWITCH OFF AND REMOVE BATTERIES BEFORE CLEANING.
- DO NOT IMMERSE IN WATER.
- DO NOT USE ABRASIVE CLEANERS OR SHARP UTENSILS TO CLEAN THIS APPLIANCE.
- WIPE SURFACES WITH A SOFT DRY CLOTH. DRY THOROUGHLY BEFORE FURTHER USE.

- STORE THE APPLIANCE IN A COOL, DRY PLACE.

BATTERY DISPOSAL

- SPENT OR EXPIRED BUTTON CELLS MUST BE DISPOSED SAFELY AWAY FROM CHILDREN ACCESS AND RECYCLED IN COMPLIANCE WITH YOUR LOCAL AUTHORITY REGULATIONS.
- FOLLOW THE BATTERY MANUFACTURER'S SAFETY, USAGE, AND DISPOSAL INSTRUCTIONS.

DISPOSAL

- Coopers of Stortford use recyclable or recycled packaging where possible.
- Please dispose of all packaging, paper, cartons in accordance with your local authority recycling regulations.
- At the end of the product's lifespan please check with your local council authorised household waste recycling centre for disposal.
- Plastics, polybags – this is made from the following recyclable plastic.

Code & Symbol	 OTHER
Type of Plastic	Please Specify (includes ABS, Nylon)

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UK: The Enterprise Department Ltd.