Rechargeable Air Compression Leg Massagers

Item No: P138

Instructions for Use v001: 30/07/24

Please read and retain these instructions for future reference

Java
wellness

ENGLISH INSTRUCTIONS

Power Details:

Rated voltage: 3.7V

Rated power: 5w

Input : DC5V 1.2A

Battery capacity: 1800 MA

Accessories:

1x Massager

1x Type-C cable

Please read this instruction manual thoroughly before starting and ensure you are familiar with the operation of your new Rechargeable Air Compression Leg Massager.

SAFETY PRECAUTIONS

Please start from the first gear strength and the first gear temperature. If you feel that the strength is too high or the temperature is too high, please change at the appropriate gear in time. Please use it in the non-charged state of the product;

Press the power key for a long time when not in use;

If you want to stop using it during use, please turn off the power first;

Do not use this product when driving;

Non-professional personnel shall not disassemble or repair without permission;

Do not use if the power cord is damaged;

Do not use when water or other liquids enter the product;

In case of any discomfort, please stop using this product immediately;

Do not place the product in a flammable, explosive, chemical or humid place;

Do not use sharp instruments to contact the product to avoid damaging the product;

This product contains a lithium battery. The battery must be removed from the appliance before it is discarded. When removing the battery, the device must be powered off.

NOTICE:

Use with caution or under the guidance of a doctor if:

-Pregnant or postpartum women

-Suffering from bone disease-Skin allergy (e.g. scald or burn)-Patients with cancer

-People with severe heart disease and abnormal blood pressure-Epilepsy patients

-The body has implanted devices (e.g. cardiac pacemakers)-People with osteoporosis

Product characteristics

The Leg Massager produced by our company has simple and fashionable product design and convenient operation. The product massages the calf muscles by air bag extrusion. Heating has a very good soothing effect on oedema, leg muscle soreness after exercise and stiff muscles. This is your ideal home masseur.

Action mechanism

The Massage Leg Massager mainly inflates and deflates the multi cavity air bag repeatedly in order to form the circulating pressure on the lower limbs and tissues, evenly and orderly squeeze the distal end to the proximal end of the lower limbs, expand blood vessels, improve local blood supply and nutrition, and improve the vitality of the tissues.

Functional features

Easy to operate and use. Operation includes temperature and mode selection.

Warm and hot compress can relieve leg muscle soreness.

Three different working modes can be adjusted independently according to the demand.

The product can be used alone or with both legs.

INSTRUCTIONS

1. Before use, please check whether the massage instrument can be started. If it cannot be started, please charge the massage instrument with the attached charging cable before use.
2. Separate the adhesive parts of the massage instrument Velcro, so that the product is fully unfolded.
3. Bind the massage instrument to the lower limbs and legs. The correct wearing direction is that the host button is upward, the digital screen is downward, and the appropriate degree of tightness is adjusted. Note: the product can be used across clothes.

4. After long pressing the power key to start, the default gear is 0.

The product enters the mode of light air pressure massage + low-grade hot compress.

The corresponding function display lights on the digital screen of the host will be on, and the green light on the screen will be on. It is timed for 15 minutes.

5. In the start-up state, press the open key gently to switch functions and modes. Press the open key to switch massage modes (leg soothing, kneading relaxation, and lifting massage). Three massage modes can be selected.

6. In the startup state, gently press the  key above to adjust the force from low to high: low gear 01 in the first gear, middle gear 02 in the second gear, and high gear 03 in the third gear. The host will shut down automatically after 15 minutes of operation.

7. In the startup state, long press 

When the heating mode is turned on, there are always three temperatures, low, medium and high, which can be adjusted according to the comfort you like. The heating mode can be used simultaneously with the massage mode.

8. After startup, the default setting time is 15 minutes. When the time is up, all loads will be turned off and enter the stop state.

WARNINGS AND SAFETY INSTRUCTIONS:

WHEN USING ELECTRICAL APPLIANCES, BASIC SAFETY PRECAUTIONS SHOULD ALWAYS BE FOLLOWED:

- READ ALL INSTRUCTIONS BEFORE USE.
- DO NOT USE OUTDOORS OR FOR COMMERCIAL USE.
- ONLY USE WITH THE MAINS ADAPTOR/CHARGER SUPPLIED.
- DO NOT USE ANY OTHER MAINS ADAPTOR/CHARGER AS THIS MAY CAUSE DAMAGE AND INVALIDATES ANY WARRANTY YOU MAY HAVE.
- BEFORE USE, CHECK THOROUGHLY FOR ANY DEFECTS AND DO NOT USE IF DEFECTS ARE FOUND. TAKE CARE NOT TO DROP THE APPLIANCE AS HEAVY IMPACTS MAY CAUSE INTERNAL DAMAGE.
- BEFORE CONNECTING THE APPLIANCE, CHECK IF THE VOLTAGE INDICATED ON THE BOTTOM OF THE APPLIANCE IS IN ACCORDANCE WITH THE MAINS VOLTAGE IN YOUR HOME.
- CHECK THE APPLIANCE AND POWER CABLE CAREFULLY FOR DAMAGE BEFORE EACH USE.
- THIS APPLIANCE CAN BE USED BY CHILDREN AGED FROM 8 YEARS AND ABOVE, AND PERSONS WITH REDUCED PHYSICAL, SENSORY OR MENTAL CAPABILITIES OR LACK OF EXPERIENCE AND KNOWLEDGE IF THEY HAVE BEEN GIVEN SUPERVISION OR INSTRUCTION CONCERNING USE OF THE APPLIANCE IN A SAFE WAY AND UNDERSTAND THE HAZARDS INVOLVED. CHILDREN SHALL NOT PLAY WITH THE APPLIANCE. CLEANING AND USER MAINTENANCE SHALL NOT BE MADE BY CHILDREN WITHOUT SUPERVISION.


- CHILDREN OF LESS THAN 3 YEARS SHOULD BE KEPT AWAY UNLESS CONTINUOUSLY SUPERVISED.
- DO NOT PLACE ON OR NEAR HEAT SOURCES.
- DO NOT COVER THE ITEM IN ANY WAY WHEN IN USE AS THIS MAY CAUSE OVERHEATING.
- DO NOT USE THIS APPLIANCE FOR ANY OTHER USE THAN THE INTENDED USE IT IS DESIGNED FOR.
- TO PROTECT AGAINST ELECTRIC SHOCK, DO NOT IMMERSE POWER CABLE, MAINS ADAPTOR/ CHARGER OR APPLIANCE, IN WATER OR ANY OTHER LIQUID.
- DO NOT OPERATE ANY APPLIANCE WITH A DAMAGED POWER CABLE OR MAINS ADAPTOR/ CHARGER OR AFTER THE APPLIANCE MALFUNCTIONS OR HAS BEEN DAMAGED IN ANY MANNER. WE RECOMMEND THAT THE POWER CABLE IS NOT KINKED OR STORED IN A TWISTED MANNER, AND THAT YOU REGULARLY CHECK ITS CONDITION.
- DO NOT LEAVE UNATTENDED.
- ANY REPAIRS MUST BE CONDUCTED BY A QUALIFIED ELECTRICIAN OR ELECTRICAL SERVICE AGENT.

CARE AND CLEANING INSTRUCTIONS

- Always switch off and unplug before cleaning.
- Do NOT immerse in water. Risk of electric shock.
- Do NOT use abrasive cleaners or sharp utensils to clean this appliance.
- Wipe surfaces with a soft dry cloth.
- Allow to dry thoroughly before further use.
- Store the appliance in a cool, dry place.

DISPOSAL

- Coopers of Stortford use recyclable or recycled packaging where possible.
- Please dispose of all packaging, paper, cartons in accordance with your local authority recycling regulations.
- At the end of the product's lifespan please check with your local council authorised household waste recycling centre for disposal.
- Plastics, polybags – this is made from the following recyclable plastic.
- At the end of the product's lifespan please check with your local council authorised household waste recycling centre for disposal.

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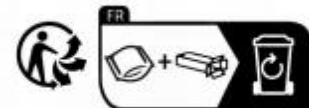Code & Symbol	
Type of Plastic	Polyethylene Terephthalate

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The Enterprise Department Ltd.