# Instructions for Use

v001: 12/09/2024

Please read and retain these instructions for future reference

## ENGLISH INSTRUCTIONS


**Power Details:**  
Input : 100-240V~  
50Hz-60Hz 0.6A MAX  
Output : 12.0V 2.0A 24.0W

Please read this instruction manual thoroughly before starting and ensure you are familiar with the operation of your new ELECTRIC FOOT & CALF MASSAGER WITH HEAT.

### Product diagram

#### 1.Product diagram
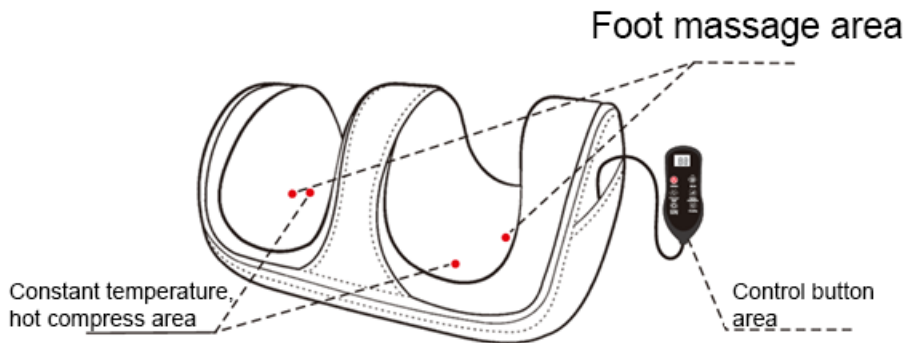
#### 2.Remote control schematic diagram

## Instructions for use

### Open key

Short press to turn on and then short press to turn off again. Turn on the red light hot compress by default, with a medium speed gear.

### Mode key

The default automatic mode for startup is to briefly press mode one, then mode two, and then mode three. Cycle in sequence.

### Forward and reverse rotation keys

After turning on the machine, press this button to exit the automatic mode, and the foot massage will switch to reverse massage. Press it again to switch to forward massage.

### Speed control key

Speed regulation: default to mid range when starting up, short press high gear, short press low gear again, and short press the third time to restore mid range.

### Heating key

When turned on, the heating function is enabled by default. Short press to turn off heating, and then short press again to turn on heating.

### Energy saving key

Energy saving design, press this button to enter energy-saving mode.

## Product maintenance and cleaning methods

### Maintenance method

- Usually, the power should be disconnected when not in use, and this product should be placed in a safe, dry, and cool place.
- Avoid contact with sharp objects to avoid piercing or damaging the appearance of the product.
- Do not store this product near high temperatures, open flames, or direct sunlight, as it may cause discoloration and damage.
- Avoid contacting this product and accessories with benzene, diluents, or any other corrosive liquids.
- When storing, do not use any power cords or wires to bundle this product or hang it with wires.

### Cleaning method

- Before cleaning, please unplug the power and gently wipe with a soft, slightly damp sponge or cloth.
- The product must not be soaked in any liquid and cannot be directly rinsed with water or any liquid.
- Do not use sharp brushes, gasoline, kerosene, furniture cleaning fluid, etc. to wipe this product.

## Common faults and solutions

Fault phenomenon	Solution
Unable to start massager	Please confirm that the power plug is correctly inserted into the socket Please confirm that the socket switch has been turned on Please confirm that the massager switch button has been pressed
Massage suddenly stops while in progress	Please check if the power plug is unplugged from the power socket or loose at the power adapter The 15 minutes set by the automatic timing device have arrived The overheat protection device is activated. Please turn off the main power supply of the product and wait for the massager to cool down

Note: If the problem persists, please do not continue to use it and contact our company's repair service center (contact information on the back cover). Do not disassemble or repair this product on your own.

### WARNINGS AND SAFETY INSTRUCTIONS:

WHEN USING ELECTRICAL APPLIANCES, BASIC SAFETY PRECAUTIONS SHOULD ALWAYS BE FOLLOWED:

- READ ALL INSTRUCTIONS BEFORE USE.
- FOR DOMESTIC INDOOR USE ONLY.
- DO NOT USE OUTDOORS OR FOR COMMERCIAL USE.
- ONLY USE WITH THE MAINS ADAPTOR/CHARGER SUPPLIED.
- **DO NOT** USE ANY OTHER MAINS ADAPTOR/CHARGER AS THIS MAY CAUSE DAMAGE AND INVALIDATES ANY WARRANTY YOU MAY HAVE.
- BEFORE USE, CHECK THOROUGHLY FOR ANY DEFECTS AND DO NOT USE IF DEFECTS ARE FOUND. TAKE CARE NOT TO DROP THE APPLIANCE AS HEAVY IMPACTS MAY CAUSE INTERNAL DAMAGE.
- BEFORE CONNECTING THE APPLIANCE, CHECK IF THE VOLTAGE INDICATED ON THE BOTTOM OF THE APPLIANCE IS IN ACCORDANCE WITH THE MAINS VOLTAGE IN YOUR HOME.
- CHECK THE APPLIANCE AND POWER CABLE CAREFULLY FOR DAMAGE BEFORE EACH USE.
- THIS APPLIANCE CAN BE USED BY CHILDREN AGED FROM 8 YEARS AND ABOVE, AND PERSONS WITH REDUCED PHYSICAL, SENSORY OR MENTAL CAPABILITIES OR LACK OF EXPERIENCE AND KNOWLEDGE IF THEY HAVE BEEN GIVEN SUPERVISION OR INSTRUCTION CONCERNING USE OF THE APPLIANCE IN A SAFE WAY AND UNDERSTAND THE HAZARDS INVOLVED. CHILDREN SHALL NOT PLAY WITH THE APPLIANCE. CLEANING AND USER MAINTENANCE SHALL NOT BE MADE BY CHILDREN WITHOUT SUPERVISION.
- CHILDREN OF LESS THAN 3 YEARS SHOULD BE KEPT AWAY UNLESS CONTINUOUSLY SUPERVISED.
- **SOME PARTS OF THIS PRODUCT CAN BECOME VERY HOT AND CAUSE BURNS. PARTICULAR ATTENTION HAS TO BE GIVEN WHERE CHILDREN AND VULNERABLE PEOPLE ARE PRESENT.**
- THE USE OF ANY ACCESSORY OR ATTACHMENTS NOT RECOMMENDED BY COOPERS OF STORTFORD MAY CAUSE INJURIES AND INVALIDATE ANY WARRANTY YOU MAY HAVE.
- USE ON A LEVEL, STABLE SURFACE.
- DO NOT COVER THE ITEM IN ANY WAY WHEN IN USE AS THIS MAY CAUSE OVERHEATING.
- DO NOT USE THIS APPLIANCE FOR ANY OTHER USE THAN THE INTENDED USE IT IS DESIGNED FOR.
- TO PROTECT AGAINST ELECTRIC SHOCK, DO NOT IMMERSER POWER CABLE, MAINS ADAPTOR/ CHARGER OR APPLIANCE, IN WATER OR ANY OTHER LIQUID.
- UNPLUG FROM OUTLET WHEN NOT IN USE AND BEFORE CLEANING.
- DO NOT OPERATE ANY APPLIANCE WITH A DAMAGED POWER CABLE OR MAINS ADAPTOR/ CHARGER OR AFTER THE APPLIANCE MALFUNCTIONS OR HAS BEEN DAMAGED IN ANY MANNER. WE RECOMMEND THAT THE POWER CABLE IS NOT KINKED OR STORED IN A TWISTED MANNER, AND THAT YOU REGULARLY CHECK ITS CONDITION.

- DO NOT LEAVE UNATTENDED.
- ALLOW TO COOL BEFORE CARRYING OR PLACING IN STORAGE.
- ANY REPAIRS MUST BE CONDUCTED BY A QUALIFIED ELECTRICIAN OR ELECTRICAL SERVICE AGENT.

### PLUG WIRING (UK & IRELAND)

This appliance is supplied with a mains adaptor/charger with BS1363 plug pins.

**DO NOT** use if the socket outlets in your home or office are not suitable for this product’s mains adaptor/ charger plug pins. Attempts to insert the mains adaptor/charger into the wrong socket is likely to cause electric hazard or damage.

**DO NOT** use with any other mains transformer/charger than the one that is supplied with this product.

### CARE AND CLEANING INSTRUCTIONS

- Always switch off and unplug before cleaning.
- Do NOT immerse in water. Risk of electric shock.
- Do NOT use abrasive cleaners or sharp utensils to clean this appliance.
- Wipe surfaces with a soft dry cloth.
- Allow to dry thoroughly before further use.
- Store the appliance in a cool, dry place.

### DISPOSAL

- Coopers of Stortford use recyclable or recycled packaging where possible.
- Please dispose of all packaging, paper, cartons in accordance with your local authority recycling regulations.
- At the end of the product’s lifespan please check with your local council authorised household waste recycling centre for disposal.
- Plastics, polybags – this is made from the following recyclable plastic.

Code & Symbol	 OTHER
Type of Plastic	Please Specify (includes ABS, Nylon)

- At the end of the product’s lifespan please check with your local council authorised household waste recycling centre for disposal.

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The Enterprise Department Ltd.