



Digital 7 Day Medicine Organiser




Thank you for purchasing our Digital 7 Day Medicine Organiser, a perfect aid to store and keep track of when you need to take your medication.

We hope you will be happy with your purchase


Please read the instruction manual carefully before use and retain for future reference.

Product Overview

1. LCD Window
2. Button Cover
3. Minute/Backlight Button
4. Set Button
5. Hour Button
6. Battery Cover
7. Strap Eyelet
8. Clip Areas

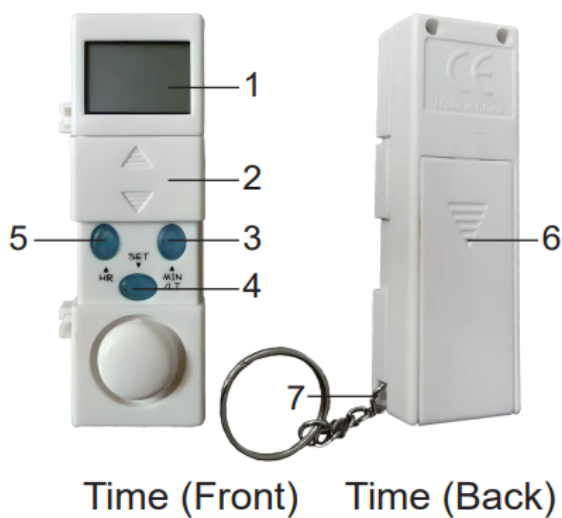


Compartment Case

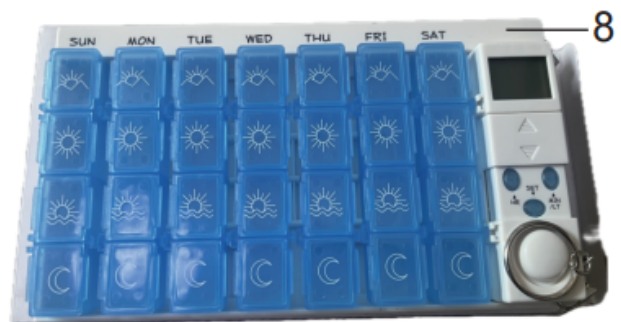

Display Information

- A. Set Pill Time Indicator
- B. Sound/Beep Indicator
- C. Clock Indicator
- D. Backlight Indicator
- E. Time

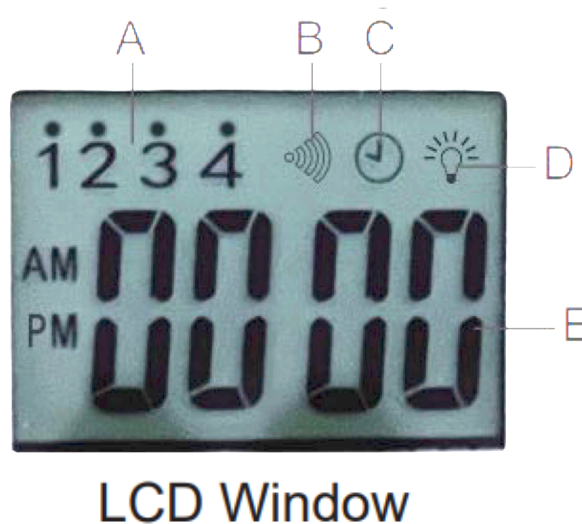

Replacing/Removing the Compartments from the Case

The pill compartments can be attached to each other or the timer via the slots/clips on the side of each item.

Battery Installation

The timer requires 2 AAA alkaline batteries (Not included).

To install, remove battery cover and insert batteries, taking note of the polarity.

Replace cover.

The timer is now ready for use. Press and hold the SET button for 3 seconds to activate timer.

Setting the Clock

Press the SET button 5 times until the time and clock indicator start to blink.

Press the Hour button until required hour is selected.

Press the Min/LT button until the required minute is selected.

Press SET button once to confirm.

To switch between 12/24 hour time, press MIN and HR at the same time.

Setting the Daily Alarm

1st time alarm (Morning)

Press the SET button once, while clock is displayed. This will show the 1st time alarm setting.

The 1 symbol will blink.

Press the HR button to select required hour, followed by the MIN/LT button to select the required minute.

Press the SET button once to confirm. A dot will now be displayed on top of the 1 symbol.

2nd time alarm (Noon)

Press the SET button twice, while clock is displayed. This will show the 2nd time alarm setting.

The 2 symbol will blink.

Press the HR button to select required hour, followed by the MIN/LT button to select the required minute.

Press the SET button once to confirm. A dot will now be displayed on top of the 2 symbol.

3rd time alarm (Evening)

Press the SET button three times, while clock is displayed. This will show the 3rd time alarm setting.

The 3 symbol will blink.

Press the HR button to select required hour, followed by the MIN/LT button to select the required minute.

Press the SET button once to confirm. A dot will now be displayed on top of the 3 symbol.

4th time alarm setting (Before Bed)

Press the SET button four times, while clock is displayed. This will show the 4th time alarm setting.

The 4 symbol will blink.

Press the HR button to select required hour, followed by the MIN/LT button to select the required minute.

Press the SET button once to confirm. A dot will now be displayed on top of the 4 symbol.

Resetting Daily Alarm

Press the SET button until the required time alarm setting is shown on the display (1, 2, 3 or 4) blinks.

Press the Hour and Min buttons simultaneously until the alarm is reset to 00:00 and the dot above the timer symbol disappears.

Press the SET button once to confirm.

LCD Backlight

Press the MIN/LT button to activate backlight. It stays on for 5 seconds.

Battery Warning:

- Do not mix alkaline, standard (carbon-zinc) and rechargeable batteries (including hydride)
- Do not mix old and new batteries
- Exhausted batteries are to be removed
- The supply terminals are not to be short-circuited
- Only batteries of the same or equivalent type as recommended are to be used
- Batteries are to be inserted with the correct polarity

Battery Disposal:

The wheellie bin symbol indicates the product and batteries must not be disposed of in the domestic waste as they can contain substances that are damaging to the environment. Please use designated collection points or recycling facilities when disposing of the item or batteries. Batteries should be removed from the product. Contact your local waste disposal agency or your retailer for further information

Product Code: P252685

Made in China for S&S Enterprises Ltd. 26 Perivale Park Middlesex UB67RL

