

SAFTTY BED RAIL WITH POCKETS

Item No: R290 / MF2087

Instructions for Use


v001: 01/07/26


COOPERS OF STORTFORD



ASSEMBLE INSTRUCTIONS:

Step by Step Instructions & Photos

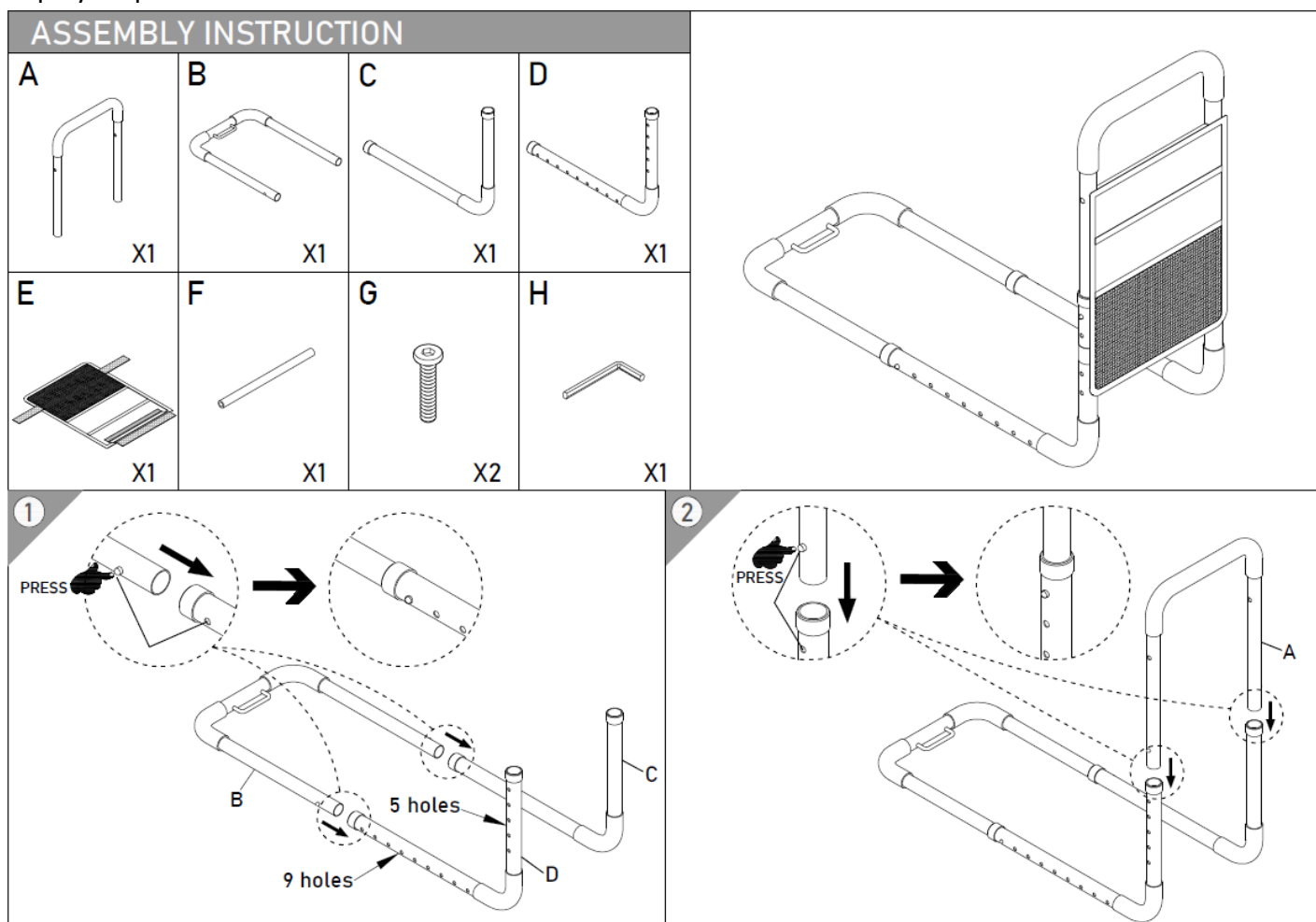

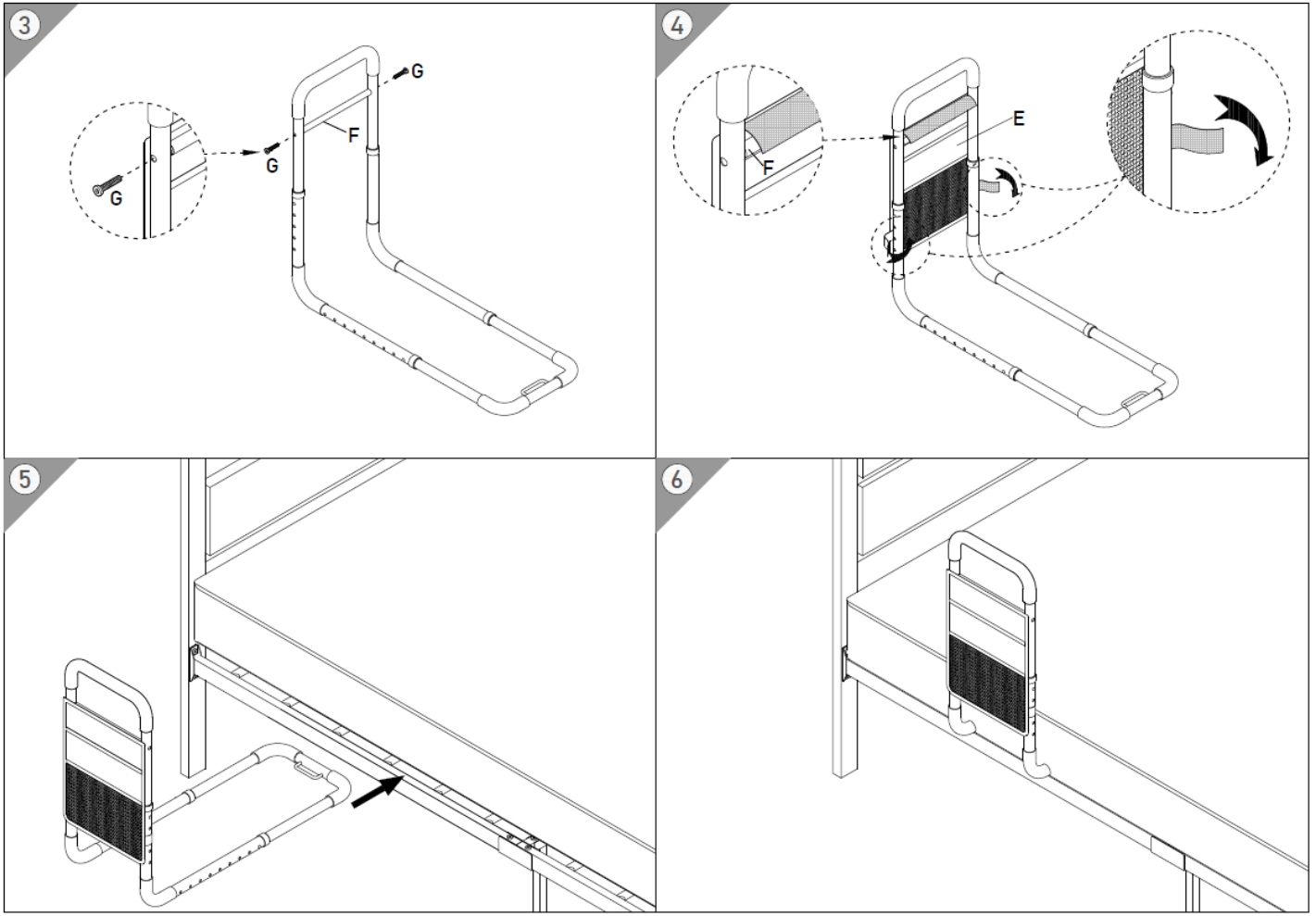

INSTRUCTIONS FOR USE:

1. Adjust the length of the bar that goes under your mattress based on the width of your mattress. The wider the mattress, the longer you adjust the bars
2. Pull up your mattress
3. Slide the product under the mattress and adjust the location where it fits your daily usage, either higher up toward your head section or lower down toward your leg section. Ideally, you can place it around your waist section
4. Adjust the height of the top bars where it is comfortable and maximize your power when using it to sit onto your bed or getting out of the bed

CARE AND CLEANING INSTRUCTIONS:

- Do NOT use abrasive cleaners or sharp utensils to clean this product.
- Store the item in a cool, dry place.
- Wipe with a clean dry cloth.

DISPOSAL

We use recyclable or recycled packaging where possible

Please dispose of all packaging in accordance with your local authority's recycling regulations

At the end of the product's lifespan check with your local council authorized household waste recycling centre for disposal.

Plastic Polybags – this is made from the following recyclable plastic

Code & Symbol	 LDPE
Type of Plastic	Low Density Polyethylene

Distributed by / Verteilt durch / Distribué par/ Distribútor / Distribuuje.

The Enterprise Department Ltd, Trading as Coopers of Stortford.

UK address: 11 Bridge Street, Bishops Stortford, Herts, UK, CM23 2JU

EU address: D.S.B. LD, 160 Bd. De Fourmies, 59100 Roubaix, France



0330 331 0300



enquiries@coopersofstortford.co.uk



<https://www.coopersofstortford.co.uk>

3Pagen

Hergestellt für 3Pagen

Sankt-Joeris-Strasse 16-28

52477 Alsdorf



+49 (0) 2404/980-204



kontakt@3pagen.de



www.3pagen.de



Points de collecte sur www.quefairedemesdechets.fr
Privilégiez la réparation ou le don de votre appareil !